

# PHAZIR™ Accessories

## Overview

The PHAZIR™ is available with an extensive range sampling attachments and options to extend the analysis capabilities and make using the PHAZIR™ even more convenient.


The PHAZIR™ Sampling adapters provide convenient method to sample liquids and solids in glass vials and quartz cuvettes. The accessories attach magnetically for quick exchange without tools. The adapters can be added to any PHAZIR™ system.

## Features

### PHAZIR™ Liquids Adapter for Vials

provides convenient sampling of liquids in 5 mm I.D. vials (1 ml) which can disposed of after use or cleaned and reused. The innovative design provides a single pass transmission measurement with 5 mm effective path length.

- Liquids sampling with disposable vials.
- 5 mm transmission path length
- Part Number: PCX-ACC1

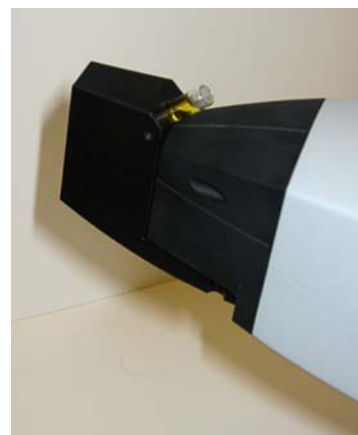

## Features

### PHAZIR™ Liquids Adapter for Cuvettes

The innovative design provides a single pass transmission measurement with 2 mm effective path length

- Liquids sampling with cuvettes
- Recommended for H<sub>2</sub>O solutions
- 5 mm transmission path length
- Part Number: PCX-ACC5



## Features

### PHAZIR™ Adapter for Diffuse Reflectance Analysis of Solids in Vials

The attachment allows diffuse reflectance measurements through the slide wall of 19 mm glass vials. The design permits vials be manually rotated in order to reposition samples and increase the effective sampling volume.

- Part Number PCX-ACC4



## Bar Code Reader

Integrated bar code reader option. Bar code values are stored within the PHAZIR.

- Part Number PCX-BarCode-1



## Spare Battery

A spare battery permits one battery to charge while the other is being used in the PHAZIR.

- Part Number: PHZR-BTRY-2



## PHAZIR™ Holster

The PHAZIR Holster provides a convenient belt pack so you can store the PHAZIR when you need your hands free. The holster has an adjustable belt which can be configured for right or left hand operation.

- Part Number: PCX-ACC2-1



## PHAZIR™ Laboratory Stand

This stand allows the PHAZIR be docked for convenient laboratory bench top operation.

- Part Number: PCX-Stand-1



## Sampling Supplies

Descriptions	Part Number	Compatibility
5 mm ID Disposable Vials, 200 vials per pack	PCX-Vial-1	PHAZIR Liquids Adapter for Disposable Vials: PN PCX-ACC1
2mm Rectangular Cuvettes, 1 vial per pack	PCX-Cuvette-2	PHAZIR Liquids Adapter for Cuvettes, PN: PCX-ACC5
15mm ID Vials, 144 vials per pack	PCX-Vial-2	PHAZIR Solids Adapter for Diffuse Reflectance measurements: PN: PCX-ACC4
Polyethylene bags, 2 x 3" , 4 mil thickness, 500 bags	PCX-Bag-1	All PHAZIR system in diffuse reflectance mode

### For additional information please contact:

Polychromix  
30 Upton Drive, Wilmington, MA 01887, USA  
Tel: +1 (978) 284-6000 Fax: +1 (978) 284-6060  
sales@Polychromix.com [www.Polychromix.com](http://www.Polychromix.com)